

Jazzy® EVO 614HD <sup>1</sup>	Group 2: Heavy Duty (K0824/K0825)
Weight Capacity	450 lbs.
Turning Radius <sup>3</sup>	22"
Overall Width <sup>3</sup>	25.4"
Overall Length <sup>3</sup>	40.5"
Maximum Speed <sup>2</sup>	Up to 4.88 mph
Ground Clearance <sup>3</sup>	3"
Front Casters	6"
Drive Wheels	13.75" perforated
Rear Casters	6"
Drivetrain	Two in-line HD motors
Braking System	Regenerative and electromechanical
Suspension Type	Full independent front and rear Active-Trac® ATX Suspension
High-back Seat Sizes	18" x 18-20" 20" x 18-20" 22" x 20-22"
Seat-to-Floor Range	22-24"
Per Charge Range <sup>2,5</sup>	Up to 13 miles at 450 lbs. Up to 21.5 miles at 200 lbs.
Standard Electronics	60A Dynamic Linx
Component Weights <sup>6</sup>	Base: 125.7 lbs. Seat: 57.8 lbs. Batteries: (2) 36 lbs. each;
Total Weight <sup>4</sup>	255.5 lbs.
Battery Requirements	(2) NF-22
Battery Charger	5 amp, off-board
Warranty	Frame: 5-year limited Drivetrain: 13-month limited Electronics: 13-month limited Battery: 13-month limited

- <sup>1</sup> Pride® FDA Class II Medical Devices are designed to aid individuals with mobility impairments.
- <sup>2</sup> Varies with user weight, terrain, battery type, battery condition and tire pressure.
- <sup>3</sup> Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or - 3%).
- <sup>4</sup> Total weight includes base, seat and batteries.
- <sup>5</sup> Tested in accordance with ANSI/RESNA, WC Vol. 2, Section 4 and ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Testing conducted at maximum weight capacities. This specification can be subject to a variance of +10%, -5%.
- <sup>6</sup> Battery weight may vary based on manufacturer.

**PRIDE**  
MOBILITY  
LIVE YOUR BEST®

(US) 800-800-8586 • (Canada) 888-570-1113

pridemobility.com |    

HCPCS codes provided should not be considered as legal advice and do not guarantee reimbursement. DME providers are responsible for determining appropriate billing codes when submitting for insurance reimbursement.

The information contained herein is correct at time of publication; we reserve the right to alter specifications without prior notice.

NOTE: All specifications subject to change without notice. The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3%. Special order seats or non-standard dimensions will not be fold-flat and depth-adjustable.

0005460-041923

**Jazzy**<sup>®</sup>  
**EVO 614HD**

**OUR MOST ADVANCED  
HEAVY DUTY POWER  
CHAIR YET!**



FDA Class II  
Medical Device<sup>1</sup>

**AMERICA'S FAVORITE POWER CHAIR**

pridemobility.com |    

# OUR MOST ADVANCED HEAVY DUTY POWER CHAIR YET!

Height, width, depth and angle adjustable armrests

Optional swing-away joystick

Depth-adjustable, limited recline, memory foam high-back seat with swivel

Removable color-through shrouds

Heavy duty seat post

13.75" perforated drive wheels

Access to controller and batteries without the user needing to exit the chair



Front and rear independent Active-Trac® ATX Suspension


FDA Class II Medical Device<sup>1</sup>

## The Future of Heavy Duty Mobility

The Jazzy® EVO 614HD has tall 13.75-inch drive wheels, the centerpiece for Pride's beloved Mid-Wheel 6® platform. Active-Trac® ATX Suspension adapts to the terrain, so grassy fields and gravel paths feel smoother than they ever have before. And now, with a new memory foam seat, your Jazzy EVO 614HD is ready to provide all-day comfort right out of the box.

- 3" ground clearance
- 13.75" perforated mid-wheel drive tires
- Front and rear independent Active-Trac® ATX Suspension
- Removable matte-finish color-through shrouds
- Depth-adjustable, limited recline, memory foam heavy duty single seat post allows the user to swivel
- Access to controller and batteries without the user needing to exit the chair
- Standard lap belt (60" - 18" or 20" W, 70" - 20" W)

## Available Colors

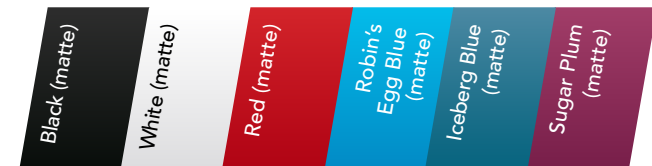

## Accessories

- Cane/crutch holder
- Cell phone holder
- Cup holder
- Flag kit
- Oxygen tank holder
- Rear basket
- Saddlebag
- Stealth® cushions
- Walker holder
- Weather cover
- 2.1A XLR USB charger
- 50-90" lap belts

## Options

- Elevated leg rests
- Solid seat pan, limited recline, high-back with swivel
- Heel loops
- Swing-away joystick
- Stealth leg rests
- Unoccupied transit
- Residual limb support



Front and rear independent Active-Trac® ATX Suspension



Large footplate with removable rubber cover for easy cleaning



High-back memory foam seat

## Military Patches

